



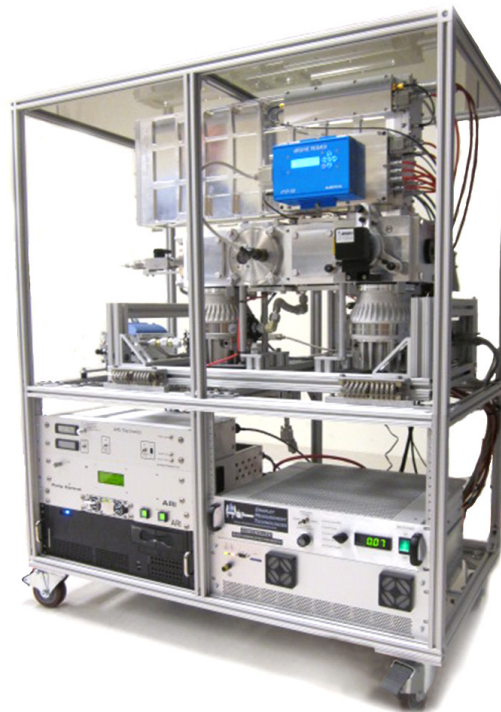

## AMS Aerosol Mass Spectrometer Systems

*Measure real-time, non-refractory, size-resolved particulate chemical composition and mass.*

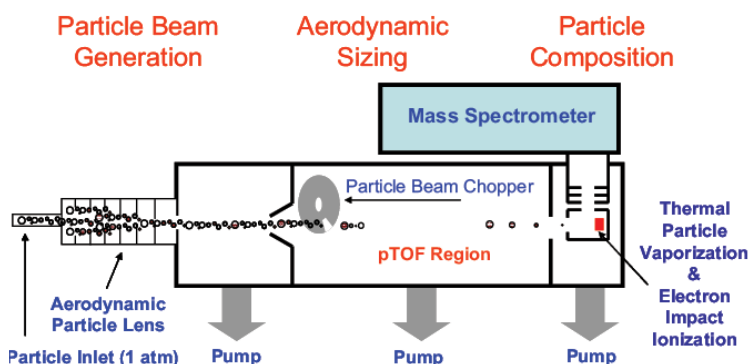

### APPLICATIONS

- Climate change and air quality research.
- Organic aerosol quantification and analysis.
  - Separation and quantification of organic components including HOA (hydrocarbon-like organic aerosol, linked to primary combustion sources) and OOA (oxygenated organic aerosol, linked to secondary aerosol sources).
  - Elemental composition (O:C, H:C).
- Mobile measurements from ship, truck and aircraft platforms.
- Fast response plume studies up to 100 Hz.
- Aerosol chamber studies.
- Combustion exhaust monitoring and source characterization.

### ADVANTAGES

- Particle beam source for efficient separation of gas and particle.
- Thermal particle vaporization with electron impact ionization source.
- Direct linear detection of sulfate, nitrate, ammonium, chloride and organic aerosol species.
- Fast response, up to 100 Hz mass spectra.
- Particle aerodynamic diameter determined from particle time-of-flight (velocity) measurements using a particle beam chopping technique.
- Several mass spectrometers to choose from: quadrupole, compact, and high resolution TOF systems.



## AMS

### SPECIFICATIONS:

Detection Limit (S/N =3) dependent on mass spectrometer option:

Mass Spectrometer System	Detection Limit* (ng/m <sup>3</sup> )	Mass Resolving Power (m/Δm)(m/z)	Mass Range
C-ToF-AMS	1.2	800	1-1000
HR-ToF-AMS:(V-mode)	2.9	2500	1-1200
:(W-mode)	32	5000	1-1200

\*Detection limits are for 1-minute integration, 3σ. Detection limits depend on chemical species. Typical values for nitrate are listed (organic DL is ~10x higher, sulfate DL is ~2x higher and ammonium DL is ~20x higher).

**Particle Size Range:** 0.04 to ~ 1.0 micrometers

**Data Rate:** 1-5 minute typical data reporting interval.  
Maximum mass spectra data rate 100 Hz (ToF MS systems only).  
Maximum size distribution data rate 150 Hz.

**Data Format:** Custom acquisition and analysis software for mass loadings and size distributions. Specialized routines for high resolution data analysis (O:C ratios).

**Sample Flow:** 0.85 l min<sup>-1</sup>

**Available Options:** Light scattering probe, black carbon detection module.

**Size/Weight:** 41" x 24" x 53", 385 lbs  
[104.14 cm x 83.82 cm x 134.62 cm, 175 kg]

**Electric Power:** 600 W; 110VAC/60Hz or 220VAC/50Hz

### REFERENCES

Jayne, J.T., D.C. Leard, X. Zhang, P. Davidovits, K.A. Smith, C.E. Kolb, and D.R. Worsnop, Development of an Aerosol Mass Spectrometer for Size and Composition. *Analysis of Submicron Particles, Aerosol Science and Technology*, 33, 49-70, 2000.

Jimenez, J.L., J.T. Jayne, Q. Shi, C.E. Kolb, D.R. Worsnop, I. Yourshaw, J.H. Seinfeld, R.C. Flagan, X. Zhang, K.A. Smith, J. Morris, and P. Davidovits, Ambient Aerosol Sampling with an Aerosol Mass Spectrometer. *Journal of Geophysical Research - Atmospheres*, 108(D7), 8425, doi:10.1029/2001JD001213, 2003.

DeCarlo, P.F., J.R. Kimmel, A. Trimborn, M.J. Northway, J.T. Jayne, A.C. Aiken, M. Gonin, K. Fuhrer, T. Horvath, K. Docherty, D.R. Worsnop, and J.L. Jimenez, Field-Deployable, High-Resolution, Time-of-Flight Aerosol Mass Spectrometer, *Analytical Chemistry*, 78: 8281-8289, 2006.

Canagaratna, M.R., J.T. Jayne, J.L. Jimenez, J.D. Allan, M.R. Alfarra, Q. Zhang, T.B. Onasch, F. Drewnick, H. Coe, A. Middlebrook, A. Delia, L.R. Williams, A.M. Trimborn, M.J. Northway, P.F. DeCarlo, C.E. Kolb, P. Davidovits, D.R. Worsnop, Chemical and Microphysical Characterization of Ambient Aerosols with the Aerodyne Aerosol Mass Spectrometer, *Mass Spectrometry Reviews*, 26, 185- 222, 2007.

Drewnick, F., S.S. Hings, P.F. DeCarlo, J.T. Jayne, M. Gonin, K. Fuhrer, S. Weimer, J.L. Jimenez, K.L. Demerjian, S. Borrmann, D.R. Worsnop. A new Time-of-Flight Aerosol Mass Spectrometer (ToF-AMS) – Instrument Description and First Field Deployment, *Aerosol Science and Technology*, 39:637- 658, 2005