




AERODYNE RESEARCH, Inc.

Dual Laser HNO₃/HONO TILDAS

Sensitive, rapid, highly specific and continuous measurements of HNO₃ and HONO in ambient air.

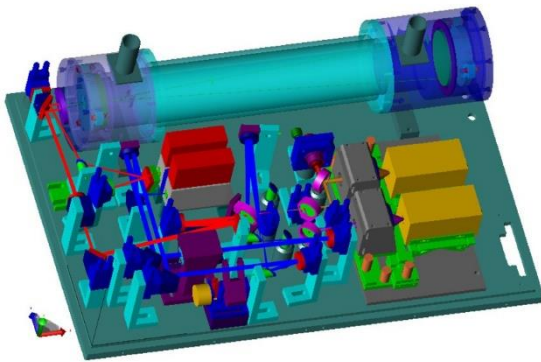

APPLICATIONS

- Extremely sensitive detection of nitric acid (HNO₃) and nitrous acid (HONO) to studies of ambient air quality, environmental nitrogen deposition/production, and others.
- Fast flow allows for <1 s time response
- Can be used with Aerodyne Inertial inlet and Active Passivation to reduce sample-wall interactions
- Eddy Covariance measurements.
- Fast response plume studies.
- Air quality monitoring.
- Mobile measurements from ship, truck, and aircraft platforms.

ADVANTAGES

- Absolute trace gas concentrations without calibration gases.
- Fast time response.
- Free from interferences by other atmospheric gases or water vapor.
- Turnkey and unattended operation.
- Ready to be deployed in field measurements and on moving platforms.
- Two lasers allow simultaneous measurement of more species.
- Optical pathlength of either 76 meters or 210 meters.

Performance Specifications (76 meter cell)

	HNO ₃	HONO
Precision at 1 sec	0.2 ppb	0.4 ppb
Precision at 100 sec	0.07 ppb	0.15 ppb


The precision is the standard deviation of a data stream measured at the typical ambient mixing ratio. Accuracy before calibration is typically 2%.

Enhanced Measurement Options

Inertial inlet for particle separation with fast time response (see right)

Multiple valve control for calibration/zeroing at inertial inlet

Active passivation to improve time response to <1 s

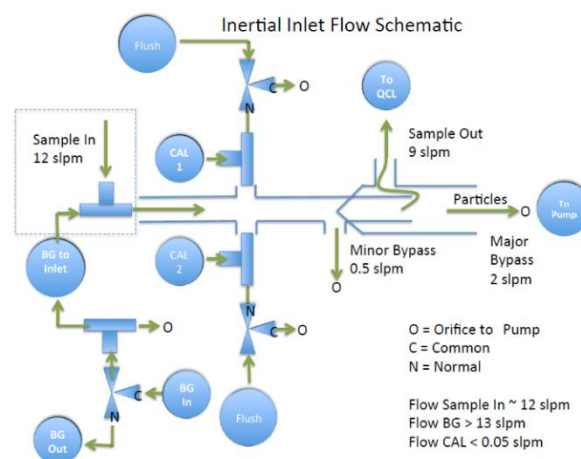

Instrument components

Core instrument
Thermoelectric chiller
Keyboard, mouse, and monitor
Vacuum pump (customer specified)

Instrument Operating Conditions

Operating temperature: 10 to 35 °C
Sample flow rate: 0 to 20 slpm

Data Outputs

RS-232, USB drive, ethernet

MECHANICAL SPECIFICATIONS FOR DUAL LASER TRACE GAS MONITOR:

Dimensions: 560 mm x 770 mm x 640 mm (WxDxH)
Weight: 75 kg
Electrical Power: 250-500 W, 120/240 V, 55/60 Hz (without pump)

MULTIPASS CELL:

Choice of 76 meter standard cell (V=0.5 liters) or 210 meter "Super Cell" (V=2liters)

REFERENCES:

Nelson, D.D. et al., Optics Let. 31, 2012-2014, 2006.
McManus, J.B. et al., Applied Physics B, DOI: 10.1007/s00340-006-2407-7 (2006).
McManus, J.B., M.S. Zahniser, D.D. Nelson, L.R. Williams, and C.E. Kolb, Infrared laser spectrometer with balanced absorption for measurements of isotopic ratios of carbon gases., Spectrochim. Acta A, 58, 2465-2479, (2002).
McManus, J.B., D.D. Nelson, J.H. Shorter, R. Jiménez, S. Herndon, S. Saleska, and M.S. Zahniser, A high precision pulsed QCL spectrometer for measurements of stable isotopes of carbon dioxide, J. Modern Optics, 52, 2309-2321 (2005).
Saleska, SR; J. Shorter, S. Herndon, R. Jimenez, B. McManus, D. Nelson, M. Zahniser, What are the instrumentation requirements for measuring the isotopic composition of net ecosystem exchange of CO₂ using eddy covariance methods? Isotopes in Environmental and Health Studies, 42 (1), 117 (2006).
Nelson, D.D., J. B. McManus, S. C. Herndon, M. S. Zahniser, B. Tuzson and L. Emmenegger, New Method for Isotopic Ratio Measurements of Atmospheric Carbon Dioxide Using a 4.3 μm Pulsed Quantum Cascade Laser, Appl. Phys. B 90, 301–309 (2008).
Tuzson, B, J. Mohn, M. J. Zeeman, R. A. Werner, W. Eugster, M. S. Zahniser, D. D. Nelson, J. B. McManus, L. Emmenegger, High precision and continuous field measurements of δ¹³C and δ¹⁸O in carbon dioxide with a cryogen-free QCLAS, Appl. Phys. B 92, 451-458 (2008).